



Weck  
**Auto Endo5**  
 Automatic Clip Applier

**Greater confidence**

The Weck Auto Endo5 Automatic Clip Applier provides greater confidence in secure occlusion using Weck Hem-o-lok Polymer Ligation Clips. Its innovative design keeps each clip seated in the applier jaws during application and delivers enhanced clip retention on the chosen structure.

**Secure-locking clip technology**

Hem-o-lok Clips deliver unparalleled ligation security compared to leading competitors,\* thanks to their integrated ridges and proprietary distal locking mechanism. The M/L clip has a clip aperture of 4.9 mm and a closed clip length of 9.2 mm.

**Procedural flexibility**

With its ergonomic design, longer shaft, and 15 Hem-o-lok M/L Clips – the only clips indicated for ligating up to a 10 mm vessel or tissue structure through a 5 mm port – the Auto Endo5 Clip Applier increases procedural flexibility.

**Benefits**



**The Clinician**

The Hem-o-lok Clip's integrated ridges and distal locking mechanism provide unparalleled clip security compared to leading competitors.\*



**Your Institution**


Cost-savings opportunity compared to leading auto endo 5 mm metal clip appliers.



**The Patient**

Flexible clip hinge enables ligation of up to a 10 mm structure through a 5 mm port, helping to minimise patient trauma.

## Weck Auto Endo5 Automatic Clip Applier Features

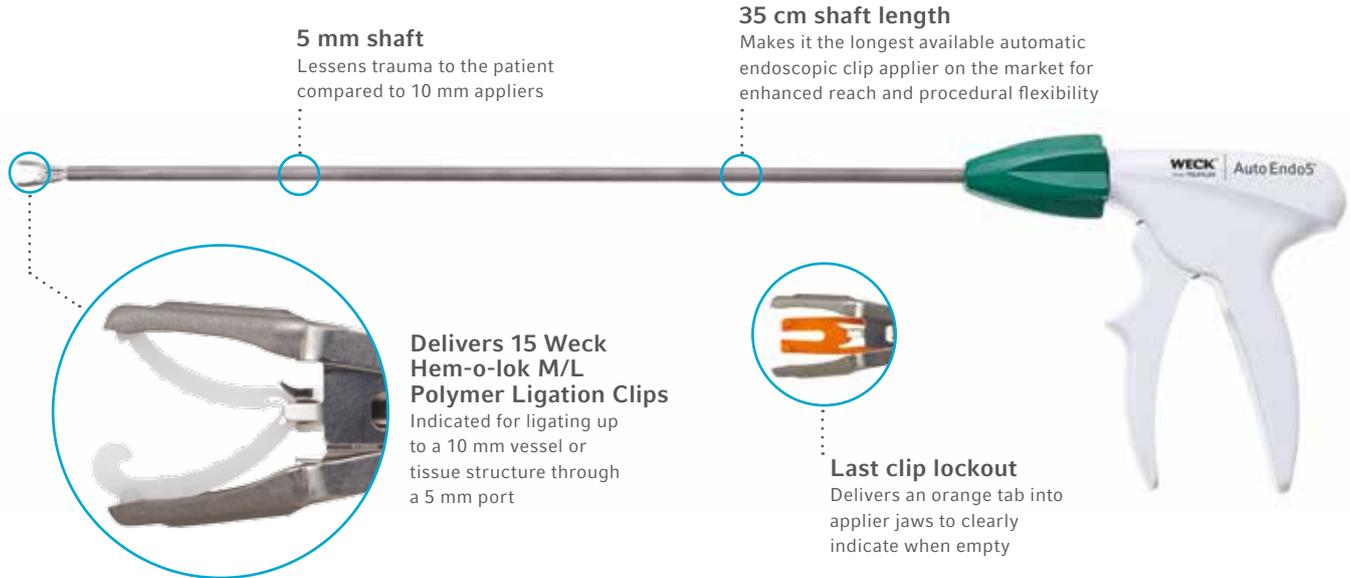

## Weck Hem-o-lok Clip Features

**Integrated ridges**  
designed to prevent slippage with a reliable 360° grip

**Locking mechanism**  
provides secure closure with tactile and auditory feedback


**Three points of contact**  
in the applier jaws for enhanced stability

**Flexible hinge**  
enables 10 mm structure ligation through a 5 mm port

**Bow shape**  
for easy removal with appropriate instrument

### Product information

PRODUCT CODE	PRODUCT DESCRIPTION	QUANTITY/SALES UNIT
AE05ML	Weck Auto Endo5 5 mm Automatic Endoscopic 35 cm Clip Applier	3
544121T	Weck Hem-o-lok Reusable Endoscopic Take a Part Clip Remover	1

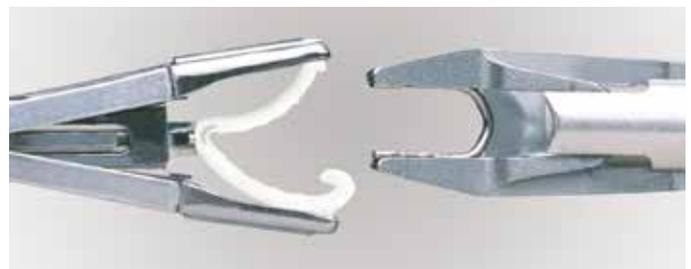

With the Hem-o-lok Clip's large distal opening and polymer flexibility, the Auto Endo5 Clip Applier is indicated for ligating larger structures than comparable 5 mm metal clip appliers.

### Applicable specialties

- General
- Gynecologic
- Thoracic
- Bariatric
- Paediatric
- Cardiac
- Colorectal
- Urologic
- ENT

**For further information, contact your Teleflex surgical sales professional.**

To place an order, contact customer service at +353 (0)9 06 46 08 00.

## Don't just clip it, lok it.

\*Data on file (2013 internal study), Teleflex Incorporated, Report #D001591. Based on ex-vivo testing of three common clip failure modes conducted on porcine carotids, sample size = 33, p ≤ 0.05. Clinical performance cannot be extrapolated from the data. Testing pressures range beyond physiological pressures. The Weck Hem-o-lok medium/large clip was compared to Ethicon LIGAMAX 5, Ethicon 10 mm LIGACLIP, and Covidien Endo Clip III. The Weck Hem-o-lok large clip was compared to Ethicon 12 mm LIGACLIP and Covidien Endo Clip. Hem-o-lok Ligating Clips are contraindicated for use as a fallopian contraceptive tubal occlusion device and contraindicated for use in ligating the renal artery during laparoscopic donor nephrectomies.

Ethicon, LIGACLIP, and LIGAMAX 5 are trademarks or registered trademarks of Johnson & Johnson. Covidien and Endo Clip are registered trademarks of Covidien. Teleflex, the Teleflex logo, Endo5, Hem-o-lok, and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this catalogue may not be available in all countries. Please contact your local representative. All data current at time of printing (09/2016). Subject to technical changes without further notice. ©2016 Teleflex Incorporated. All rights reserved. 94 30 54-00 00 01 REV A 09 16 02

**Teleflex**

Distributed by:  
Teleflex Headquarters International, Ireland · Teleflex Medical Europe Ltd. · IDA Business & Technology Park  
Dublin Road · Athlone · Co Westmeath · Tel. +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.int@teleflex.com  
United Kingdom Tel. +44 (0)14 94 53 27 61 · info.uk@teleflex.com  
South Africa Tel. +27 (0)11 807 4887 · assist.africa@teleflex.com